

Code No. 27362

Human Mac-2 binding protein (Mac-2bp) Assay Kit - IBL
INTRODUCTION

Mac-2 binding protein (Mac-2bp), known as 90K, is a highly N-glycosylated, secreted protein, identified as a ligand of Galectin-3. It is considered that through interaction with Galectin-3, Mac-2bp promotes homotypic cell-cell contact or regulates cell adhesion. And it has been reported that Mac-2bp levels in blood have associations with various human cancers or several viral infectious diseases. This ELISA kit can measure concentration of Mac-2bp.

PRINCIPLE

This kit is a solid phase sandwich ELISA using 2 kinds of highly specific antibodies. Tetra Methyl Benzidine (TMB) is used as a coloring agent (Chromogen). The strength of coloring is proportional to the amount of Human Mac-2bp.

MEASUREMENT RANGE

0.78 - 100 ng/mL

INTENDED USE
For research use only, not for use in diagnostic procedures.

- This IBL's assay kit is capable for assay of human Mac-2bp in serum, EDTA-plasma and cell culture media.
- The guide line of dilution rate for serum and plasma samples is from 500 to 1,000-fold.

KIT COMPONENT

1	Precoated plate	:	Anti-Human Mac-2bp (8A2) Mouse IgG MoAbAffinity Purify	96Well x 1
2	Labeled antibody Conc.	:	(30X) HRP conjugated Anti-Human Mac-2bp (67A1) Mouse IgG MoAb Affinity Purify	0.4mL x 1
3	Standard	:	Recombinant Human Mac-2 binding protein	0.5mL x 2
4	EIA buffer	:	1% BSA, 0.05% Tween20 in PBS	50mL x 1
5	Solution for Labeled antibody	:	1% BSA, 0.05% Tween20 in PBS	12mL x 1
6	Chromogen	:	TMB solution	15mL x 1
7	Stop solution	:	1N H ₂ SO ₄	12mL x 1
8	Wash buffer Conc.	:	(40X) Phosphate buffer	50mL x 1

OPERATION MANUAL
1. Materials needed but not supplied

- Plate reader (450nm)
- Graduated cylinder and beaker
- Refrigerator (as 4°C)
- Paper towel
- Washing bottle for precoated plate
- Disposable test tube for "2, Labeled antibody Conc." and "6, Chromogen"
- Micropipette and tip
- Deionized water
- Graph (semilogarithmic) paper
- Tube for dilution of Standard


2. Preparation

- 1) Preparation of wash buffer
"8, Wash buffer Conc." is a concentrated (40X) buffer. Adjust the temperature of "8, Washing buffer Conc." to room temperature and then, mix it gently and completely before use. Dilute 50 mL of "8, Wash buffer Conc." with 1,950 mL of deionized water and mix it. This is the wash buffer for use. This prepared wash buffer shall be stored in refrigerator and used within 2 weeks after dilution.
- 2) Preparation of Labeled antibody
"2, Labeled antibody Conc." is a concentrated (30X). Dilute "2, Labeled antibody Conc." with "5, Solution for Labeled antibody" in 30 times according to required quantity into a disposable test tube. Use this resulting solution as Labeled antibody.
Example)
In case you use one strip (8 well), the required quantity of Labeled antibody is 800 µL. (Dilute 30 µL of "2, Labeled antibody Conc." with 870 µL of "5, Solution for Labeled antibody" and mix it. And use the resulting solution by 100 µL in each well.)
This operation should be done just before applying labeled antibody.
The remaining "2, Labeled antibody Conc." should be stored at 4°C in firmly sealed vial.
- 3) Preparation of Standard
Put just 0.5 mL of deionized water into the vial of "3, Standard" and mix it gently and completely. This solution is 200 ng/mL human Mac-2bp standard.
- 4) Dilution of Standard
Prepare 8 tubes for dilution of "3, Standard". Put 230 µL each of "4, EIA buffer" into the tube.
Specify the following concentration of each tube."

Tube-1	100 ng/mL
Tube-2	50 ng/mL
Tube-3	25 ng/mL
Tube-4	12.5 ng/mL
Tube-5	6.25 ng/mL
Tube-6	3.13 ng/mL
Tube-7	1.56 ng/mL
Tube-8	0.78 ng/mL
Tube-9	0 ng/mL (Test Sample Blank)

Put 230 µL of Standard solution into tube-1 and mix it gently. Then, put 230 µL of tube-1 mixture into tube-2. Dilute two times standard solution in series to set up 8 points of diluted standard between 100 ng/mL and 0.78 ng/mL. Tube-9 is the test sample blank as 0 ng/mL.

See following picture.

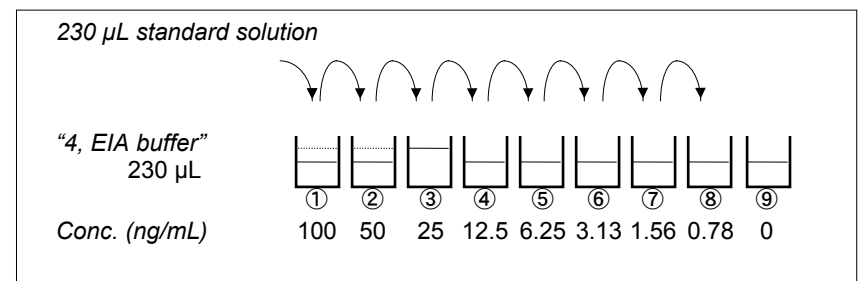
5) Dilution of test sample

Test samples should be diluted with "4, EIA buffer" suitably.

Serum or plasma samples have to be diluted with "4, EIA buffer" accordingly. The recommended dilution for them is from 500 to 1,000-fold. In case of the absorbance of sample is over than the assay range, it is necessary to dilute it more.

<Example of 500-fold dilution of serum or plasma>

1. Add 20 µL of serum or plasma to 380 µL of "4, EIA buffer" in a tube and mix them well.
2. Pipette 20 µL of 20-fold diluted serum or plasma from the tube of above first dilution and add it to 480 µL of "4, EIA buffer" in another tube, and mix them well.
(In the case of 1,000-fold dilution, pipette 10 µL of 20-fold diluted sample and add it to 490 µL of "4, EIA buffer")
3. This 500-fold diluted serum or plasma should be applied as a test sample according to the measurement procedure.

3. Measurement procedure

All reagents shall be brought to room temperature approximately 30 minutes before use. Then mix it gently and completely before use. Make sure of no change in quality of the reagents. Standard curve shall be prepared simultaneously with the measurement of test samples.

Reagents	Test Sample	Standard	Test Sample Blank	Reagent Blank
	Test sample 100 µL	Diluted standard (Tube 1-8) 100 µL	EIA buffer (Tube-9) 100 µL	EIA buffer 100 µL
Incubation for 60 minutes at 4 °C with plate lid				
Washing 4 times				
Labeled Antibody	100 µL	100 µL	100 µL	-
Incubation for 30 minutes at 4 °C with plate lid				
Washing 5 times				
Chromogen	100 µL	100 µL	100 µL	100 µL
Incubation for 30 minutes at room temperature (shielded)				
Stop solution	100 µL	100 µL	100 µL	100 µL
Read the plate at 450nm against a Reagent Blank within 30 minutes after addition of Stop solution.				

- 1) Determine wells for reagent blank. Put 100 µL each of "4, EIA buffer" into the wells.
- 2) Determine wells for test sample blank, test sample and diluted standard. Then, put 100 µL each of test sample blank (tube-9), test sample and dilutions of standard (tube-1-8) into the appropriate wells.
- 3) Incubate the precoated plate for 60 minutes at 4 °C after covering it with plate lid.
- 4) Wash each well of the precoated plate 4 times with wash buffer using a washing bottle or a plate washer in following way.
After shaking off (or aspiration of) the solution in wells, fill each well with wash buffer and shake off the wash buffer completely from the precoated plate. This procedure must be repeated 4 times. Then, drain the precoated plate completely on paper towel.
In case of using a plate washer, we recommend manually washing in the manner mentioned above at the last one time.
Please refer to 5) and 6) in SPECIAL ATTENTION below, and be careful not to miss a well.
- 5) Pipette 100 µL of labeled antibody solution into the wells of test samples, diluted standard and test sample blank.
- 6) Incubate the precoated plate for 30 minutes at 4°C after covering it with plate lid.
- 7) Wash the precoated plate 5 times in the same manner as 4).
- 8) Take the required quantity of "6, Chromogen" and put it into a disposable test tube. Then, pipette 100 µL from the test tube into every well. Please do not return the rest of used chromogen in the test tube into "6, Chromogen" bottle in order to avoid contamination.
- 9) Incubate the precoated plate for 30 minutes at room temperature in the dark. The solution of Chromogen will turn blue.
- 10) Add 100 µL of "7, Stop solution" to all wells. Mix the solution by tapping the side of precoated plate. The solution will turn yellow by addition of "7, Stop solution".
- 11) Remove any dirt or drop of water on the bottom of the precoated plate and confirm there is no bubble on the surface of the solution. Then, run the plate reader and conduct measurement at 450 nm against a reagent blank. The measurement shall be done within 30 minutes after addition of "7, Stop solution".

SPECIAL ATTENTION

- 1) Test samples should be measured soon after collection. For the storage of test samples, store them frozen and do not repeat freeze/thaw cycles. Thaw

the test samples at a low temperature and mix them completely before measurement.


- 2) Test samples should be diluted with "4, EIA buffer", suitably.
- 3) Duplicate measurement of test samples and standard is recommended.
- 4) Use test samples in neutral pH range. The contaminations of organic solvent may affect the measurement.
- 5) Use only wash buffer in this kit for washing the precoated plate. Insufficient washing may lead to the failure in measurement.
- 6) Remove the wash buffer completely by tapping the precoated plate on paper towel. Do not wipe wells with paper towel.
- 7) "6, Chromogen" should be stored in the dark due to its sensitivity against light. Avoid contact of Chromogen with metals.
- 8) Measurement should be done within 30 minutes after addition of "7, Stop solution".

CALCULATION OF TEST RESULT

Subtract the absorbance of test sample blank from all data, including standards and unknown samples before plotting. On a semilogarithmic paper the concentration of the standards (x-axis, logarithmic) are plotted against their corresponding absorbance (y-axis, linear). Draw the best smooth curve through these points. Read the concentration for unknown samples from the standard curve. In automated method, 4 parameter logistics can generally give a good fit.

Example of standard curve




Conc. (ng/mL)	Absorbance (450nm)
100	1.929
50	1.691
25	1.363
12.5	0.983
6.25	0.664
3.13	0.399
1.56	0.236
0.78	0.133
0 (Test Sample Blank)	0.004

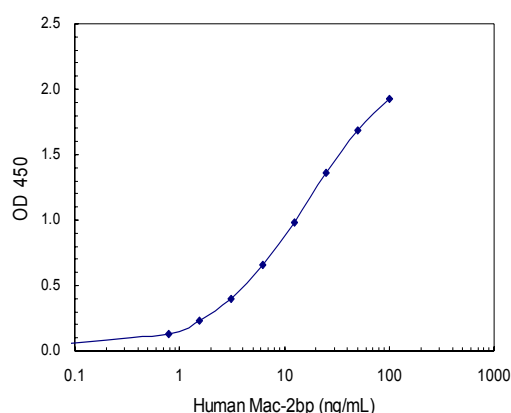

The typical standard curve is shown above. This curve can not be used to derive test results. Please run a standard curve for each assay.

PERFORMANCE CHARACTERISTICS

1. Dilution linearity

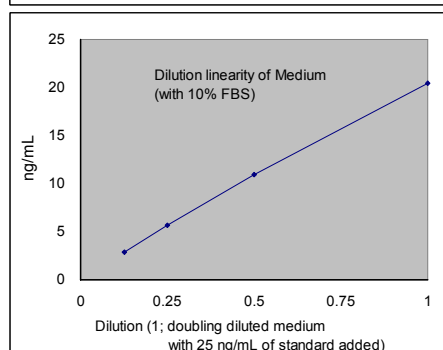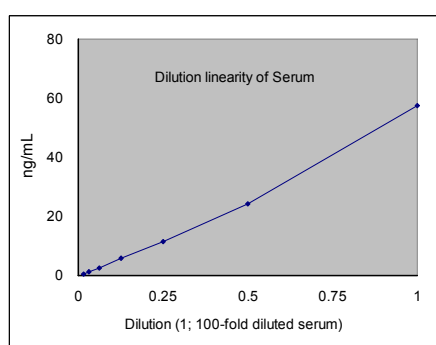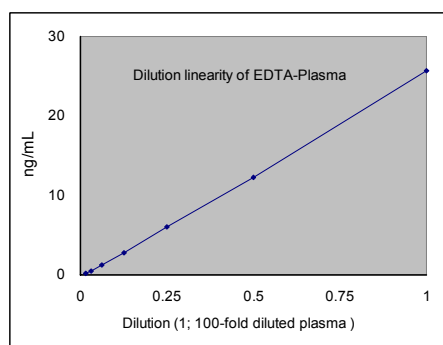

2. Added Recovery Assay

Specimen	Additive Amount (ng/mL)	Theoretical Value (ng/mL)	Measured Value (ng/mL)	%
Human Plasma (EDTA) (x500)	12.5	17.52	15.98	91.2
	6.25	11.27	10.34	91.7
	3.13	8.15	7.28	89.4
Human Serum (x500)	12.5	17.11	16.11	94.1
	6.25	10.86	9.71	89.3
	3.13	7.74	6.90	89.2
Medium with 10% FBS (x10)	12.5	12.5	11.16	89.3
	6.25	6.25	5.33	85.3
	3.13	3.13	2.77	88.5

3. Intra - Assay

Mean Value (ng/mL)	SD (ng/mL)	CV (%)	n
42.25	2.46	5.8	22
11.83	0.53	4.5	22
2.96	0.10	3.4	22

4. Inter - Assay

Mean Value (ng/mL)	SD (ng/mL)	CV (%)	n
39.38	3.21	8.2	7
10.88	0.66	6.1	7
2.85	0.11	3.8	7

5. Sensitivity

0.08 ng/mL

The sensitivity for this kit was determined using the guidelines under the National Committee for Clinical Laboratory Standards (NCCLS) Evaluation Protocols. (National Committee for Clinical Laboratory Standards Evaluation Protocols, SC1, (1989) Villanova, PA: NCCLS.)

PRECAUTION FOR INTENDED USE AND/OR HANDLING

1. All reagents should be stored at 2 - 8°C. All reagents shall be brought to room temperature approximately 30 minutes before use.
2. "3, Standard" is lyophilized products. Be careful to open this vial.
3. "7, Stop solution" is a strong acid substance. Therefore, be careful not to have your skin and clothes contact "7, Stop solution" and pay attention to the disposal of "7, Stop solution".
4. Dispose used materials after rinsing them with large quantity of water.
5. Precipitation may occur in "2, Labeled antibody Conc.", "4, EIA buffer" or "8, Wash buffer Conc.", however, there is no problem in the performance.
6. Wash hands after handling reagents.
7. Do not mix the reagents with the reagents from a different lot or kit.
8. Do not use expired reagents.
9. This kit is for research purpose only. Do not use for clinical diagnosis.

STORAGE AND THE TERM OF VALIDITY

Storage Condition : 2 - 8°C

The expiry date is specified on outer box.

REFERENCE

1. Cheung KJ, Tilleman K, Deforce D, Colle I, Van Vlierberghe H. The HCV serum proteome: a search for fibrosis protein markers. *J Viral Hepat.* 2009 Jun;16(6):418-29.
2. Ohshima S, Kuchen S, Seemayer CA, Kyburz D, Hirt A, Klinzing S, Michel BA, Gay RE, Liu FT, Gay S, Neidhart M. Galectin 3 and its binding protein in rheumatoid arthritis. *Arthritis Rheum.* 2003 Oct;48(10):2788-95.
3. Marchetti A, Tinari N, Buttitta F, Chella A, Angeletti CA, Sacco R, Mucilli F, Ullrich A, Iacobelli S. Expression of 90K (Mac-2 BP) correlates with distant metastasis and predicts survival in stage I non-small cell lung cancer patients. *Cancer Res.* 2002 May 1;62(9):2535-9.
4. Matarrese P, Fusco O, Tinari N, Natoli C, Liu FT, Semeraro ML, Malorni W, Iacobelli S. Galectin-3 overexpression protects from apoptosis by improving cell adhesion properties. *Int J Cancer.* 2000 Feb 15;85(4):545-54.
5. Inohara H, Akahani S, Kohts K, Raz A. Interactions between galectin-3 and Mac-2-binding protein mediate cell-cell adhesion. *Cancer Res.* 1996 Oct 1;56(19):4530-4.

Version 1.

Made in Japan.